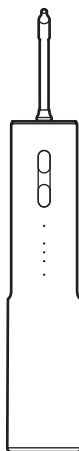
Panasonic

Operating Instructions (Household use) Oral Irrigator

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Panasonic

English

Operating Instructions (Household use) Oral Irrigator Model No. EW-NJ8A

Contents

Safety precautions	4
Intended use.....	11
Parts identification	13
Charging	15
How to use.....	18
Tips of use	22
Cleaning	23
Optional parts (Sold separately).....	26
Troubleshooting.....	27
Battery life.....	30
Removing and disposing of the built-in rechargeable battery.....	30
Specifications	32

Thank you for purchasing this Panasonic product.
**Before operating this unit, please read these instructions
completely and save them for future use.**

The USB power adaptor is sold separately.

Safety precautions

Make sure to follow these instructions.

To prevent accidents, injuries or damage to property, please follow the instructions below.

- The following chart indicates the degree of damage caused by improper operation.

	DANGER	Indicates a potential hazard that will result in serious injury or death.
	WARNING	Indicates a potential hazard that could result in serious injury or death.
	CAUTION	Indicates a potential hazard that could result in minor injury or property damage.

WARNING

► Power supply

- **Fully insert the power plug of the USB power adaptor securely.**

- Failure to do so may cause fire or electric shock.

- **Do not use in a way exceeding the rating of the household outlet or the wiring.**

- Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

► About the USB power adaptor and USB cable

- **Do not immerse them in water or wash with water.**

- **Never use the appliance if they are damaged or if the power plug fits loosely in a household outlet.**

- **Do not damage them.**

- Do not deface, modify, forcefully bend or pull them. Also, do not place anything heavy on, twist or pinch them.

- Doing so may cause electric shock or fire due to a short circuit.

- **Do not connect or disconnect them with a wet hand.**

- Doing so may cause electric shock or injury.

- **Disconnect the USB plug when cleaning the main unit.**

- Failure to do so may cause electric shock or injury.

- **Regularly clean them to prevent dust from accumulating.**

- Failure to do so may cause fire due to insulation failure caused by humidity.

Remove the USB power adaptor and wipe it with a dry cloth.

WARNING

- **Do not charge this product using any USB cable other than the supplied USB cable. Also, do not use any other product with the supplied USB cable.**
 - Doing so may cause fire, electric shock, or injury.
 - The USB cable provides the overheat protection function. If charging is abnormal or the USB plug (Type-C) heats up, disconnect the power supply, check the socket [A ⑨], and remove any water or foreign objects from the socket.
- **Always plug the USB cable to a power supply with a rated voltage of 5.0 V \equiv and a rated output current of 1.0 A which meets requirements of the Safety Extra Low Voltage (SELV) system.**

- Failure to do so may cause burn or fire due to a short circuit.

It is necessary to separately purchase a power adaptor if required to power this product. The Panasonic power adaptor ES-RU1-02 is recommended.

► In case of an abnormality or malfunction

- **Immediately stop using and remove the adaptor if there is an abnormality or malfunction.**
 - Failure to do so may cause fire, electric shock, or injury.

<Abnormality or malfunction cases>

- **The main body, USB power adaptor or USB cable is deformed or abnormally hot.**
- **The main body, USB power adaptor or USB cable smells of burning.**
- **There is abnormal sound during use or charging of the main body.**
 - Immediately request inspection or repair at an authorized service center.

WARNING

English

► This product

- **This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, or leave in a high temperature environment.**
 - Doing so may cause overheating, ignition, or explosion.
- **Do not modify or repair.**
 - Doing so may cause fire, electric shock, or injury.
 - Contact an authorized service center for repair (battery change, etc.).
- **Never disassemble except when disposing of the product.**
 - Doing so may cause fire, electric shock, or injury.

► Preventing accidents

- **Do not store within the reach of children or infants. Do not let them use it.**
 - Doing so may cause an accident or injury due to accidental ingestion of removable parts.
- **This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.**
- **Children should be supervised to ensure that they do not play with the appliance.**
 - Failure to do so may cause an accident or injury.

CAUTION

- ▶ **Preventing damage to teeth and gums**
- **People who are unable to operate should not use this appliance.**
- **People with weak oral senses should not use this appliance.**
 - Doing so may cause injury or damage to teeth and gums.
- **People who may have periodontal disease, have teeth being treated, or concerned about symptoms within their mouth should consult a dentist before use.**
 - Failure to do so may cause injury or damage to teeth and gums.
- ▶ **This product**
- **Never use for a purpose other than oral cleaning.**
 - Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, tongue, or throat.
- **Do not fill the tank with water warmer than 40 °C.**
 - Doing so may cause burns.
- **Do not share nozzle with your family or other people.**
 - Doing so may cause infection or inflammation.
- ▶ **Note the following precautions**
- **Do not drop or subject to shock.**
 - Doing so may cause injury.
- **Do not allow metal objects or trash to adhere to the power plug of the USB power adaptor or the plugs of the USB cable.**
 - Doing so may cause electric shock or fire due to a short circuit.

CAUTION

English

- **Do not charge when the appliance socket is wet.**
 - Doing so may cause electric shock or fire.
- **Do not wrap the USB cable around the power adaptor before storage.**
 - Doing so may cause the wire in the USB cable to break with the load, and may cause fire due to a short circuit.
- **Disconnect the USB power adaptor from the household outlet when not charging.**
 - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.
- **After use, clean appliance and the nozzle and dry them.**
 - Failure to do so may cause damage to your health due to propagation of mold and bacteria.

Handling of the removed battery when disposing

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- **Never peel off the tube.**
 - Doing so may cause overheating, ignition, or explosion.

WARNING

- **After removing the rechargeable battery, do not keep it within reach of children and infants.**
 - The battery will harm the body if accidentally swallowed. If this happens, consult a physician immediately.
- **If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.**
 - Failure to do so may cause eye injury.
 - Consult a physician immediately.

CAUTION

- **If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.**
 - Failure to do so may cause inflammation.

Intended use

- Please use this product together with a toothbrush. It is hard to remove tartar adhering to the tooth surface only by a water jet, so it is necessary to brush your teeth with a toothbrush. We recommend that you rinse your mouth with a water jet after brushing your teeth.
- When you first begin using the oral irrigator, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the oral irrigator for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the oral irrigator and consult your dentist.
- Panasonic's specified adaptor with splash-resistant design and a USB cable with an overheating protection feature should be used to ensure safety. The adaptor is not completely waterproof, so prevent it from getting wet.
- Do not put salt, mouthwash or other chemicals into the water tank. Doing so may cause malfunction.
- Do not dry the main unit with a dryer or a fan heater. Doing so may cause malfunction or deformation of the parts.
- Oral irrigators shall be used only with tap water or lukewarm water, or specific solutions as defined in the manufacturer's instructions.

- Children shall not play with the appliance.
- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.
- Charge the appliance correctly according to these Operating Instructions. (See page 15 “Charging”).
- This appliance contains battery that is only replaceable by skilled persons. Please contact an authorized service center for the details of repair.

This symbol on the product means that not generated higher voltage than the safety extra low voltage (SELV). In which protection against electric shock relies on power supply.

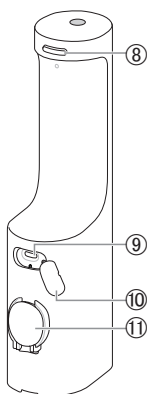
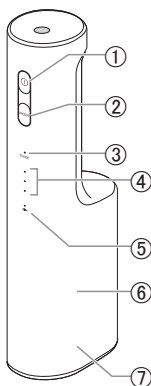
Parts identification

English


A Front

Back


(While the socket cover is open)


B

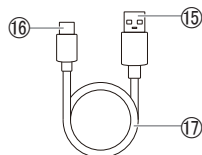C

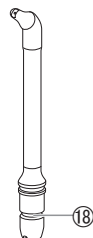D

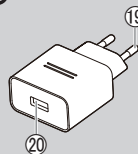Tips:

There may be moisture inside the main body (e.g. in the water tank), left from the distilled water used for product testing. This does not pose any hygiene concerns.

A Main unit

- ① Power switch (OFF/ON)
- ② MODE (Mode/Level) selector button
- ③ GUIDE (Mode indicator) (blue)
- ④ Water pressure level indicator (blue)
- ⑤  Charge status indicator (red)
- ⑥ Water tank (Inside)
- ⑦ Main body
- ⑧ Nozzle release button
- ⑨ Socket (Type-C)
- ⑩ Socket cover
- ⑪ Water tank cap
- ⑫ Cover plate
- ⑬ Bottom cover
- ⑭ Bottom cover release button

B USB cable

- ⑮ USB plug (Type-A)
- ⑯ USB plug (Type-C)
- ⑰ Cord

C Nano bubble nozzle plus

- The nozzle is a consumable.
- ⑱ Nozzle joint

D USB power adaptor for general use (ES-RU1-02)
(Sold seperately)

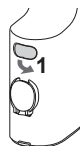
- ⑲ Power plug
- ⑳ USB port (Type-A)

Charging


English

- Please wipe off any water droplets from the socket or socket cover.


- 1 Open the socket cover on the back of the main body.

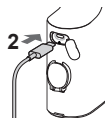- 2 Plug the Type-C plug of the USB cable into the socket on the main body.

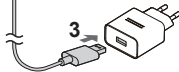
- 3 Insert the Type-A plug into a power adaptor with a rated voltage of 5.0 V --- and a rated output current of 1.0 A which meets requirements of the Safety Extra Low Voltage (SELV) system.


- Check the direction of the Type-A plug terminals, hold the Type-A plug, and insert or pull it out straight. (Deflection or insertion in a wrong direction may cause deformation or malfunction.)

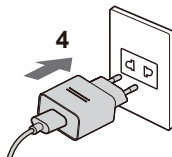- 4 Insert the USB power adaptor into a household outlet.

5 After charging, unplug the USB cable or/and the power adaptor from the household outlet to ensure safety and save energy.

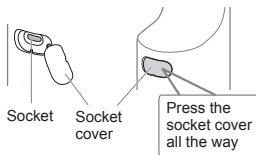- Remove the Type-C plug from the socket, and close the socket cover.
- It takes about 3 hours to fully charge this product. (The charging time may differ depending on the remaining battery level.)
- Even if the charging time exceeds 3 hours, the product will not malfunction.
- Keep the socket cover closed when the product is not charged, especially while you are using the product or washing the main body with water.
- Be sure to close the socket cover when using or washing.

■ Operating time after a full charge

Approx. 20 minutes*

(When the continuous water jet mode 3 is set)

* It may vary the ambient temperature and usage conditions.

■ About the charging reminder function

If the battery is running low,  (the charge status indicator) flashes (once/sec) as a reminder under the following circumstances:

- When the power switch is turned on
- When the MODE selector button is pressed
- When the power switch is turned off (the indicator will light off after 5 seconds)

→ Please charge the product when the charging reminder appears.

Indicator display during charging

		
Charging in progress	Charging complete	Charging abnormal (Flashing twice a second)

(● : Lights on / ○ : Lights off) (✳ : Flashing)

How to check if charging is completed

After charging is completed, if you remove and then insert the USB cable back,  (the charge status indicator) lights on and then lights off in about 5 seconds.

► Notes

- The recommended ambient temperature for charging is 10 °C - 35 °C. It may take longer to charge, or  will flash rapidly and the battery may fail to charge properly in extremely high or low temperatures.
- The charging time may differ when the appliance is used for the first time, or if it has not been used for more than 6 months.
In addition,  (the charge status indicator) may not light on in the first several minutes of charging, but it will light on as the charging proceeds.
- The battery will deteriorate if not used for over 6 months, so make sure to perform a full charge at least once every 6 months.

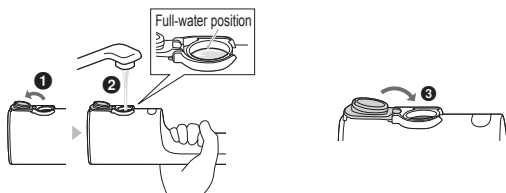
How to use

- Check that the socket cover is closed.


1 Fill the water tank with water.

- ① Open the water supply cover.
- ② Add water while keeping the main unit level.
- ③ Close the water supply cover securely.


■ If the water supply cover comes off

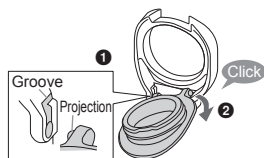- ① Align the projection on one side of the water tank cap with the groove on the main body, and then press it all the way down until you hear a click.
- ② Insert the other side in the same way.

2 Install the nozzle.

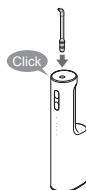- Fully install the nozzle until you hear a click.

3 Press the MODE selector button, and select your desired mode/level.

- About the mode and level:
Each time you press the MODE selector button, the mode is switched once in the following order:

Continuous water jet mode

1 → 2 → 3 → Interdental

guide mode 1 → 2 → 3 →

Continuous water jet mode 1...

- About the two water jet modes:

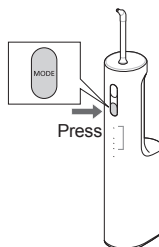
Continuous water jet mode and interdental guide mode.

The continuous water jet mode is designed to clean all areas of the oral cavity with a steady water jet.

The interdental guide mode is designed to alternate between strong and weak water jets. When the water jet is strong, direct the water jet onto the space between the teeth; when the water jet is weak, move the water jet to the next interdental area.

- About the GUIDE (mode indicator):

In the continuous water jet mode, the GUIDE (mode indicator) lights off; in the interdental guide mode, the GUIDE (mode indicator) lights on.

Mode/Level	Continuous water jet mode			Interdent guide mode		
	1	2	3	1	2	3
GUIDE (Mode indicator)	○	○	○	●	●	●
Water pressure level indicator (Strong)  (Weak)	○	○	●	○	○	●
	○	●	●	○	●	●
	●	●	●	●	●	●

(Lights on ●
Lights off ○)

■ About the mode memory function

The selected mode is remembered. From the second use, the product shows the mode that was selected last time.

- When the product has not been used for a long time or when the battery is used up, the mode will reset to the continuous water jet mode 1.

4 Clean the inside of your mouth with a water jet.

- 1 Place the nozzle in your mouth.
 - 2 Close your mouth gently.
 - 3 Turn on the power switch.
- During use, you can press the MODE selector button to switch modes.

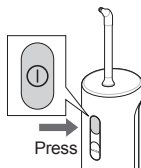
■ Operation time with full water

Approx. 80 seconds.

(When the continuous water jet mode 3 is set and the product is used in an upright position)

- If the main body is tilted during use, the operation time may be reduced.
- If the main unit is tilted too much, the water does not come out from the nozzle.

5 To prevent the water from splashing about, turn off the power switch while keeping the nozzle in your mouth.

► Notes

- The appliance does not work during charging.
- Be sure to use clean tap water or lukewarm water. (Do not put salt, mouthwash, medication, or other similar substances into the water tank.)
- If the water feels cold, please use lukewarm water with a temperature of 40 °C or less.
- When using the appliance for the first time, please press the power switch and check the water pressure level while water is coming out. (When using the product for the first time or if the sensation feels too intense, please start with the continuous water jet mode 1 or the interdental guide mode 1.)
- The appliance may not operate in an ambient temperature approximately 5 °C or lower.
- Do not push the tip of the nozzle too hard against the teeth or gums.
- Do not turn on the power switch when the water tank is empty unless you are drying the moisture inside the main unit. (Doing so may cause malfunction.)

Tips of use

Adjust the angle and movement of the nozzle according to the purpose and the area you want to clean.

Point

For areas where the gums are sensitive, please use the weak water pressure level 1.

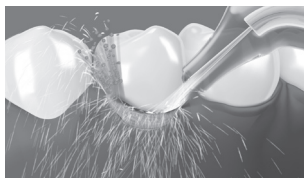
► Interdental cleaning

- Direct the water jet onto the space between the teeth.
- Direct from the front and rear sides of the teeth.

► Periodontal pocket cleaning

- Aim the water jet at the space between teeth and gums.
- Slowly shift the water jet along the gums.
- Clean the front and back of the teeth and behind the molars as well.

► Gums care

- Stimulate the gums with the water jet.

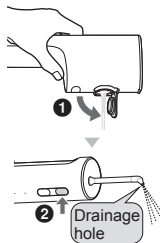
Cleaning

- After use, drain off any residual water in the water tank.
- Clean the product after each use.

1 Pour out the water from the water tank, and drain off any remaining water from the main body.

- ① Open the water tank cap and pour out the water.
 - Be sure to drain off the water after each use.
- ② Turn on the power switch “①” to remove the remaining water.
- ③ Keep it on until no water comes out.

2 Remove the nozzle and bottom cover.

- ① Remove the nozzle

Nozzle release button

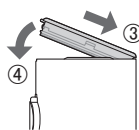
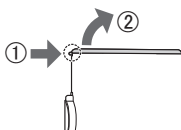
Pull out the nozzle while pressing the button.

- ② Remove the bottom cover


→ How to attach/remove the bottom cover

- To remove: ① Press and hold the bottom cover release button (○), and ② lift it up to remove it.
- To attach: Attach it in the order of ③ and ④

3 Wash the main body (as a whole), the nozzle and the inside of the water tank.

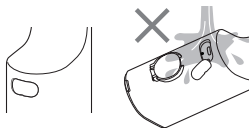- Make sure the socket cover is closed before washing with water.

Washing with water

- Avoid hitting it against the basin or other surfaces when pouring out the water.

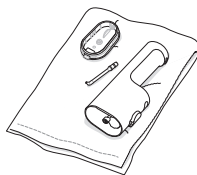
4 Store the product after drying it.

- 1 Wipe any water away from the nozzle and the main unit with a towel or the like. (Otherwise, malfunction or contamination may be caused.)


- Do not charge the product when it is wet.

- 2 Attach the bottom cover after drying it in a cool and ventilated place.

■ If the rubber ring of the bottom cover comes off

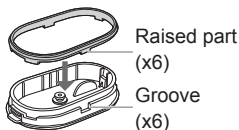With the flat side of the rubber ring facing inward, align the raised part of the rubber ring with the groove on the bottom cover, and attach the rubber ring to the bottom cover.


► Notes

- Do not clean the product with thinner, benzine, alcohol, etc. (Doing so may cause malfunction, or cracking or discolouration of the parts.)
- Use neutral kitchen detergent and rinse well with water to ensure no detergent remains after cleaning.
- Do not wash it with hot water exceeding 50 °C.

■ When the water jet is poor

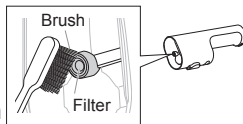Remove the bottom cover and clean any debris from the filter.
Gently scrub it clean using a soft brush or the like.


■ When it feels dirty

Remove the cover plate from the bottom cover and wash it with water.


→ How to attach/remove the cover plate

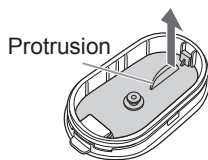- To remove: Lift it upward by pinching the raised part

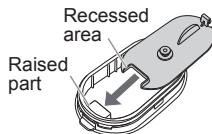- To attach: Align the raised part of the bottom cover with the recessed area of the cover plate, and then press the cover plate down until you hear a click.


■ USB cable



Before cleaning, unplug the USB cable from the power adaptor.
Wipe with a dry cloth.
(Do not rinse it with water.)

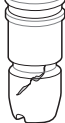Optional parts (Sold separately)

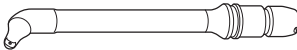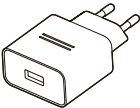- The nozzle is a consumable.
- Replacement cycle: Approx. 6 months (for hygiene concerns).
- Once the tip of the nozzle gets deformed or cracked, replace the nozzle with a new one even within 6 months.

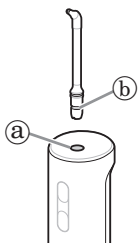
A new nozzle (tip)	A deformed (worn) nozzle (tip)	A cracked nozzle (insertion part)
		

Part number	Replacement Parts	
WEW0870	Nano bubble nozzle plus	
ES-RU1-02	USB power adaptor	

Troubleshooting

English

Problem	Possible causes and corrective measures
<div>While charging</div> <div>  <p>The charge status indicator does not light on.</p> </div> <div>  <p>flashes rapidly (twice a second).</p> </div>	<p>The appliance was just purchased or it has not been used for a while. →Keep charging. The charge status indicator will light on after a while. (See page 17.)</p> <p>The USB cable provides the overheat protection function. It automatically powers off when the Type-C plug overheats. →Disconnect the power supply, check the socket, and remove the water or foreign objects if any.</p> <p>When the socket is wet or foreign objects adhered to the socket, the USB cable will activate the self-protection function. →Disconnect the power supply, check the socket, and remove any water or foreign objects from the socket.</p> <p>The USB cable is not firmly connected to the USB power adaptor. →Check to see that the USB cable is in full contact with the USB power adaptor, and that the charge status indicator is on. (See page 15.)</p> <p>The ambient temperature is too low or too high during charging. →Please charge within the recommended temperature range of 10 °C to 35 °C.</p>

Problem		Possible causes and corrective measures
While charging	The appliance can be used for just several minutes even after a full charge.	<p>The USB cable is not firmly connected to the USB power adaptor. → Make sure that the USB plug (Type-C) of the USB cable is fully inserted into the socket of the main unit or the USB plug (Type-A) of the USB cable is fully inserted into the power supply and the charging status indicator is on. (See page 15.)</p> <p>The battery has reached the end of its life. (Approx. 3 years) → Replace it with a new one. Please contact your dealer or Service Center.</p>
	The water jet is poor.	<p>The nozzle gets clogged, deformed, or cracked. → Replace the nozzle with a new one. (See page 26.)</p> <p>The filter is clogged. → Remove the dirt from the filter. (See page 25.)</p>
While using	No water comes out.	<p>The main unit is tilted too much. → Please use the appliance with the main unit upright. (See page 20.)</p> 
	The nozzle is hard to install.	<p>It may take some time to suck up water. → Put water into the nozzle port ①, install the nozzle, and turn on the power switch.</p> 
		→ Moisten the nozzle joint ②.

Problem		Possible causes and corrective measures	
While using	The nozzle cannot be fully inserted into the main body.	▶	<p>There may be residual fragments of a damaged nozzle.</p> <p>→ While pressing the nozzle release button (C), remove the fragments with a cotton swab or the like.</p> 
	Water leaks from the bottom cover.	▶	<p>The rubber ring on the bottom cover is lifted or misaligned.</p> <p>→ Please install it properly. (See page 24.)</p>
	 does not flash before turning off.	▶	<p>The battery runs out. (Approx. 3 years)</p> <p>→ Replace it with a new one. Please contact your dealer or Service Center. (See page 30.)</p>
	It does not operate even after charging.	▶	<p>The ambient temperature is too low during use.</p> <p>→ Do not operate this product under a temperature of 5 °C or lower.</p>

If the problems still cannot be solved, contact the store where you purchased the unit or a service center authorized by Panasonic for repair.

Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge.

(The operating life of the battery may differ depending on usage or storage conditions.)

Removing and disposing of the built-in rechargeable battery

This product has a built-in rechargeable battery.

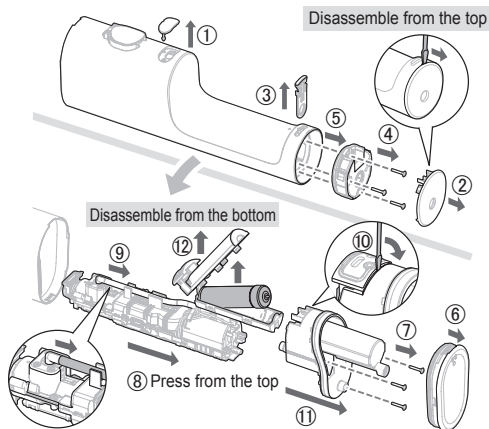
Be sure to remove the battery before disposing of the product.

- The battery is to be disposed of safely.
- This figure must only be used when disposing of the appliance, and must not be used to repair it.
- If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.
- We recommend seeking the advice of a dealer or Panasonic's authorized service center before deciding to dispose of the product, so as to extend the service life of the product as much as possible.
- If you decide to dispose of the product but do not want to disassemble it by yourself to remove the battery, we recommend taking the product to Panasonic's authorized service center which will remove and recycle the battery.

Removing the built-in rechargeable battery

- Power off the product before removing the battery. Turn on the power switch and keep the power on until the battery is completely discharged.
- Follow Steps ① ~ ⑫ as illustrated to disassemble the product and remove the battery with a screwdriver.

- During disassembly, take care not to short circuit the positive and negative terminals of the removed battery and insulate the terminals by applying tape to them.
- If you cannot remove the battery, please contact an authorized service center.

For environmental protection and recycling of materials

This appliance contains a Lithium-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Output rating of applicable adaptor (sold separately)	5.0 V \equiv 1.0 A In compliance with the Safety Extra Low Voltage (SELV) system
Charging time	Approx. 3 hours (Room temperature: 10 °C~35 °C)
Operating time	Approx. 20 minutes (at 20 °C when fully charged) (When the continuous water jet mode 3 is set)
Tank volume	Approx. 145 mL (When the continuous water jet mode 3 is set) (approx. 80 seconds)
Dimensions	Approx. 233 (H) x 53 (W) x 89 (D) mm (Not including nozzle)
Mass	Approx. 435 g (Including nozzle)

This product is intended for household use only.

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© Panasonic Corporation 2025
Panasonic Corporation