Panasonic®

Operating Instructions

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(Household) Massage Lounger

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Model No. EP-MA05
型號

English	EN2
繁體中文	TC2

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

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Contents

Before using the product

- Safety Precautions EN4
- Part names and functions of main unit EN10
- Settings..... EN11

How to use

- Before each massage (check the unit)..... EN13
- Select menu and start massage EN15
- "Auto" massage course EN16
- "Manual" massage course EN17
- Set up according to your favorite action ... EN18
- After completing the massage EN19

Maintenance

- Cleaning and maintenance EN20

Q&A

- Q&A EN21
- Troubleshooting EN22
- Sounds and sensations of the unit EN23
- Specifications EN24

The illustrations used in the instructions may be slightly different from the actual unit, and the pictures are for reference only.

Safety Precautions Always follow these instructions

This section contains instructions that you must follow to prevent personal injury or damage to property.

- The instructions are classified according to the level of injury or damage that may occur due to the mistaken use of this unit.

Please ensure that you read all instructions before using the Household Massage Lounger.

 WARNING	Denotes a potential hazard that could result in serious injury or death.
 CAUTION	Denotes a hazard that could result in minor injury or property damage.

- Instructions that you must follow are indicated using the following symbols.
(The following symbols are examples.)

 Denotes a specific operating procedure that must not be performed.	 Denotes a specific operating procedure that must be performed.
--	--

WARNING

Symptoms

- The following persons should not use this unit:
 - (1) Persons who are prohibited from receiving massages by a physician (e.g.: due to thrombosis, deep vein thrombosis of the lower limbs, pulmonary embolism, severe aneurism, acute varicose veins, phlebitis, any type of dermatitis or skin infection [including inflammation of the hypodermis], etc.)
 - (2) Persons suffering from osteoporosis, spine fracture, sprain or acute pain such as a pulled or torn muscle
 - (3) Persons who have a pacemaker, defibrillator, or other personal medical devices prone to electronic interference
 - (4) Pregnant women and women who have just given birth
 - (5) Persons with symptoms of acute low back pain, herniated disk, spondylolisthesis, spinal stenosis, degenerative lumbar spondylosis
 - (6) Persons who have abnormalities or curvature of the spine
(Otherwise it may worsen symptoms.)
- The following persons should not use the leg massage or lower back/hip massage.
 - (1) Those with symptoms of acute inflammation around the hip area, acute hip pain, piriformis muscle syndrome or sciatic neuralgia
 - (2) Those suffering from osteoarthritis of the hip and bad hip joint
 - (3) Those who feel pain, numbness, lethargy in the hip and legs during use
(Otherwise it may worsen symptoms.)
- This unit is not intended for use by persons (including children) with reduced physical sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the unit by a person responsible for their safety. Children should be supervised to ensure that they do not play with the unit.
(Otherwise it may lead to accident or injury.)

WARNING

- Persons with any of the conditions listed below or who are undergoing medical treatment should consult with their physicians before using the unit.

- (1) Persons who have a malignant tumor
- (2) Persons suffering from heart disease
- (3) Persons who have a loss of heat sensation
- (4) Persons who have sensory abnormalities due to a severe peripheral circulatory disturbance such as diabetes mellitus
- (5) Persons undergoing anticoagulation therapy
- (6) Persons who have a wound where the massager operates
- (7) Persons whose body temperature is over 38 °C (who have a fever)
(E.g.: When a person has acute inflammation symptoms (fatigued, chills, blood pressure fluctuations, etc.) or when debilitated)
- (8) Persons who require bed rest or who are in poor physical shape
- (9) Persons other than those listed above who feel unwell

- If you begin to feel unwell during use or if you do not feel effects of the massage lounger, stop using the unit immediately and consult a physician.

(Otherwise it may lead to accident, injury or become unwell.)

Take note of the following points as well

- Absolutely do not modify. Furthermore, do not attempt to disassemble or repair the unit yourself.
(Otherwise it may lead to fire, or cause the unit to function abnormally resulting in injury.)

Safety Precautions (continued)

WARNING

When in use

- When massaging the neck area, be careful of the movement of the massage rollers and avoid massaging the throat area and excessively strong massage action.
- Do not sit on the backrest or armrests.
- This unit should not be used by children. Also, do not allow children to play on the unit. Do not allow children to climb on the leg massage section, seat, backrest or armrests.
- Do not use the massage lounger while holding children.
- Do not push your hands, elbows, arms or feet against the backrest fabric covering the moving parts of the massage rollers during use.
(Otherwise it may lead to accident or injury.)

- Before use, be sure to check whether the covering fabric is damaged. (If it is damaged, please stop using it immediately and entrust repair.) Please do not use the unit when the fabric has been ripped, otherwise it may cause injury or electric shock.
- Please control the time of each massage within 15 minutes. Do not massage any one part of your body for more than 5 minutes at a time.
(Otherwise it may lead to adverse effect or injury.)
- Make sure there are no persons or pets around the unit (front, back, underneath or sides of the unit) before using it for massaging or moving the backrest or legrest. (Otherwise it may lead to accident or injury.)
- After each massage, push the power switch on the unit backside to the "O" (off) position.
(To prevent breakage or injury from misuse by children.)

If a malfunction or breakdown occurs

- If a malfunction or breakdown occurs, stop using the chair immediately and unplug the power plug from the wall socket.
(Otherwise it may lead to smoke/fire or electric shock.)

<Examples of malfunctions/breakdowns>

- Does not start when pressing .
- The power supply goes on and off when the power cord is moved.
- Burning smells or abnormal sounds occurring during operation.
- The unit becomes deformed or unusually hot.

→ Please contact customer service center for inspection and repair immediately.

WARNING

Power plug, power cord etc.

- Do not do anything that may damage the power cord or power plug.

Do not scratch, modify, use near a heater, excessively bend, twist, pull, place heavy objects on top of, or bundle the cord or plug.

(Use damaged power cord or power plug may cause electric shock, short circuit or fire.)

→ Contact the store where you bought the unit or the customer service center for repairs to the power cord or power plug.

- Do not insert or unplug the power plug with wet hands. (Otherwise it may lead to electric shock.)
- Always use sockets and wiring devices at the correct rated value.
- Use 220 V 50 Hz alternating current. (Cannot be used overseas, nor with a transformer.)
- Always insert the power plug fully into the socket. (Otherwise it may lead to electric shock or fire due to overheating.)
- Dust the power plug on a regular basis.
 - Remove dust using a dry cloth. (Otherwise it may lead to fire, due to damp shorting the insulation.)

CAUTION

Symptoms

- Persons who are otherwise healthy but have any of the conditions listed below should consult with their physicians before using the unit.

- (1) Persons whose muscles have deteriorated due to age or weight loss
 - (2) Persons with symptoms of lower back and hip pain caused by the bones and muscles, organs and nerves
 - (3) Persons who suffer bruises or sprains easily
 - (4) Persons who suffer from severe motion sickness
 - (5) Persons who have undergone heart or internal surgery in the past (Otherwise it may lead to adversely affected health.)
- If symptoms such as rashes, reddening, itching, etc. occur as a result of using this unit, stop using the unit and consult a physician. (Otherwise it may lead to accident or adversely affected health.)

Installation and movement

- The unit should not be used on top of heating appliances, such as electric carpets, etc. (Otherwise it may lead to fire.)
- Do not drag or push the unit in an installed state.
- Do not move the unit using the castors on the flooring. (Otherwise it may damage flooring.)
- Do not move with a person on it. (Otherwise it may lead to falls, causing an accident or injury.)
- Do not use in damp or humid places such as a bathroom. (Otherwise it may lead to electric shock.)

Safety Precautions (continued)

CAUTION

-  Use only on a horizontal surface.
(Otherwise it may lead to a falling accident.)
- When unplugging the power plug, hold the plug itself and do not pull on the power cord.
(Otherwise it may lead to electric shock, fire or ignition due to short circuit.)
- Before moving the unit, the angle of the leg massage section must be adjusted back to the original position.
(Otherwise it may lead to injury.)
- Place a mat or other such covering on floor when moving the unit using the castors.
(Otherwise it may damage flooring.)
- When moving the unit on its castors, remove any obstacles in the surroundings of the unit, lift up the backrest to the fully upright position, lift the legrest to a position below the waist (less than 80 cm from the floor), and move the unit slowly. (EN11)
(Otherwise the unit may fall and cause injury.)

Power supply

-  Always unplug the power plug from the wall socket when cleaning the unit or removing the cover.
(Otherwise it may lead to electric shock or burn.)
- Always unplug the power plug from the wall socket when not using.
(Since dust and humidity deteriorates insulation, which may result in fire due to electrical fault.)

Before and during use

-  Do not use the massage rollers on your head, stomach or bare skin. Also, do not place your hands or feet between the massage rollers.
- Do not clamp your knees between the leg massage section.
- Do not fall asleep while using the unit. Do not use the unit after drinking alcohol.
(Otherwise it may lead to accident or injury.)
- Do not clamp your hands or legs between the unit and the leg massage section.
- Do not clamp hands, arms or heads between the gaps of the leg massage section (legrest). For the sake of safety, please take back the leg massage section after massage.
- Do not unplug the power plug or turn the power switch "O" (off) during massage.
- Do not stand or sit on the leg massage section or the legrest. Otherwise, the end of the legrest may rotate or the unit may tilt which causes accidents or injuries.
- Do not use at the same time as another medical device.
(Otherwise it may lead to accident or deteriorated health.)




 **CAUTION**

- Check that the massage rollers are in the retracted position.
- If you want to leave the seat during massage, you must stop the massage before leaving the seat.
(Otherwise it may lead to accident or injury.)
- When retracting the leg massage section, please push it completely into position by hand.
(Otherwise it may lead to accident or injury.)
- Check that there is no foreign matter between the parts of the unit before sitting down.
 - Make sure there is no foreign matter between backrest, seat, leg massage section and pillow.
(Otherwise, hands, feet or other objects may become trapped, resulting in an accident or injury.)
- Avoid wearing anything hard on your head like a hair accessory, etc. when using the unit.
- Do not massage with any hard objects in your trouser pockets.
(Otherwise it may lead to injury.)

Take note of the following points as well

-  ● Do not spill water on the unit or controller.
(Otherwise it may lead to electric shock, fire or ignition due to short circuit.)
-  ● In the event of a power outage, unplug the power plug immediately.
(Otherwise it may lead to accident or injury when the power supply returns.)
- If you have not used the unit for a while, carefully read the operating instructions again and check that the unit is operating normally before use.
(Otherwise it may lead to malfunction or injury.)

Part names and functions of main unit

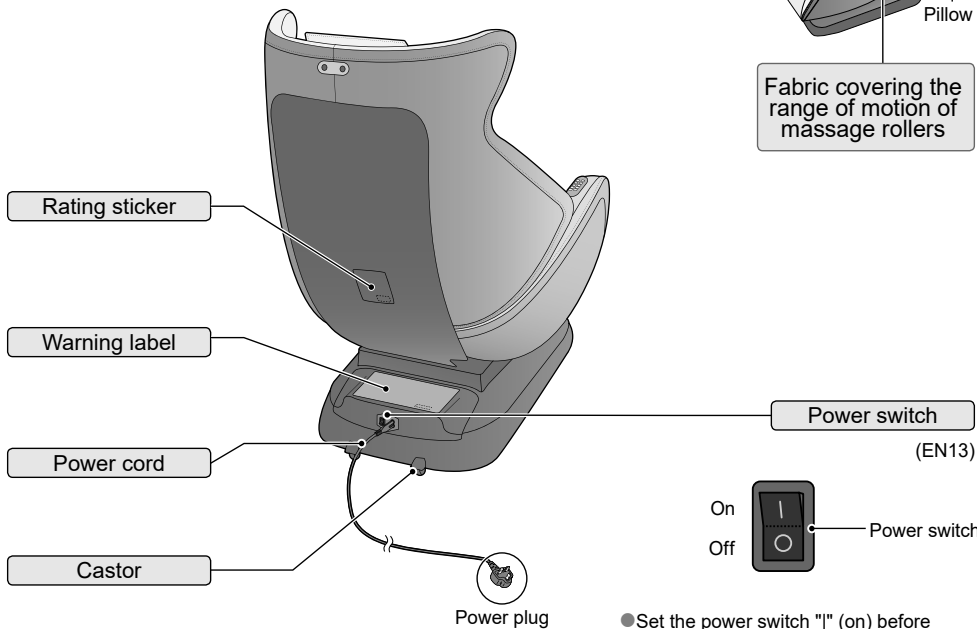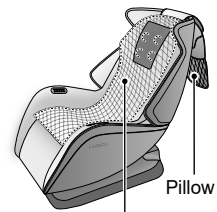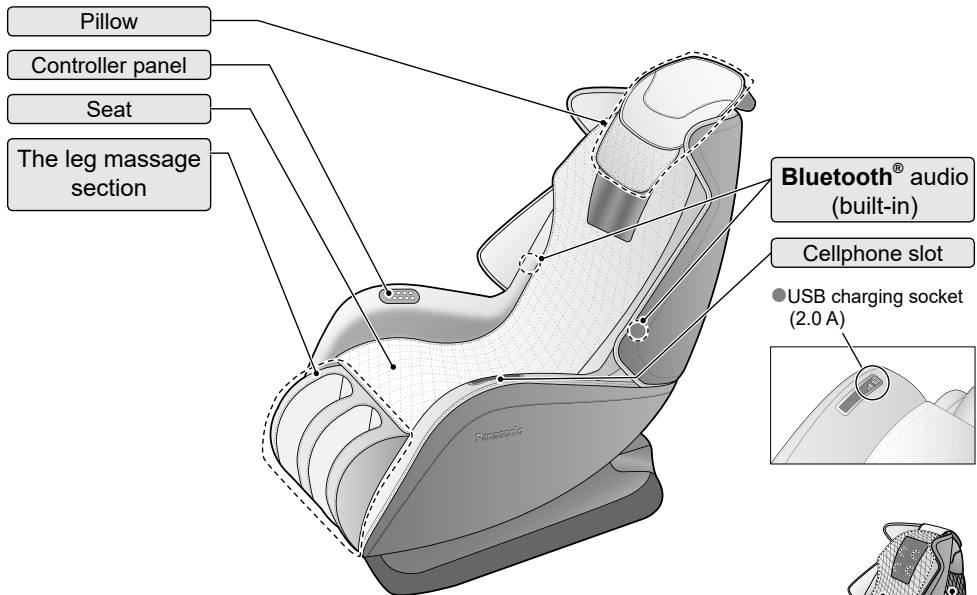EN10


Settings

- 1) Open the packing case and take out the massage lounger and accessories

Check the unit and accessories

The unit

The state of being taken out of the packing case.

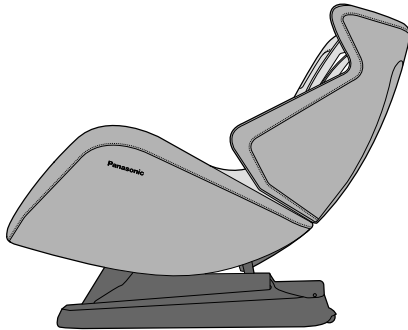Annex

- Power cord


2) Set up the unit

- The unit is heavy, so please be careful not to hurt the back and other parts during carrying. (The unit must be carried by more than 2 persons.)
The weight of the unit is approx. 43 kg.
- When you lower the unit, please take care of your feet and put it down slowly.
- Do not let go of the leg massage section until the unit body is safely laid down.

When lifting the unit and moving


If not handled correctly, the unit body will fall and damage the floor. Please fully tilt the backrest, retract the leg massage section, and then lift the massage lounger.

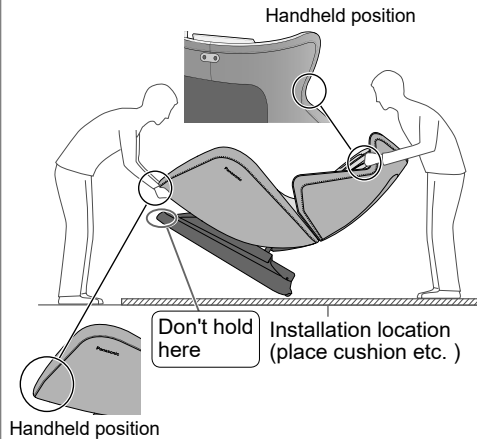Settings (continued)

3) Pay attention to the distance between the backrest and the wall when placing the unit

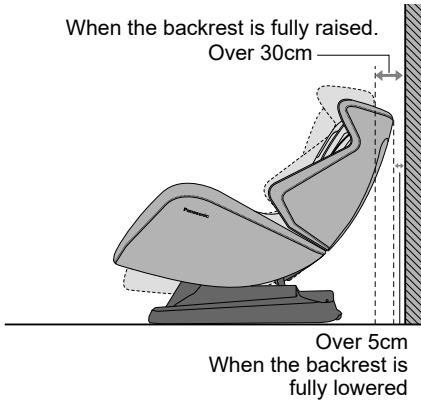
Please check whether the backrest will hit the wall

Castors can be used when moving the unit.


1 Initial use

- ① After procedures ①, ④ and ⑤ of EN13, the indicator light will flash, which means the unit enters the standby mode. The backrest of the massage lounger is automatically lifted from the delivered state to the upright state.
- ② When pressing , the indicator light is always on, which means it is at the boot mode, and at the same time, the backrest of the massage lounger tilts to the set angle (not the lowest position).
- ③ When pressing , the indicator light flashes, that is, it returns to standby mode, and the backrest of massage lounger will also return to upright state.
- ④ Turn the power switch on the back of the unit to "O" and unplug the power plug.

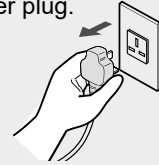2 Moving


Lift the unit to a position no higher than the position of waist

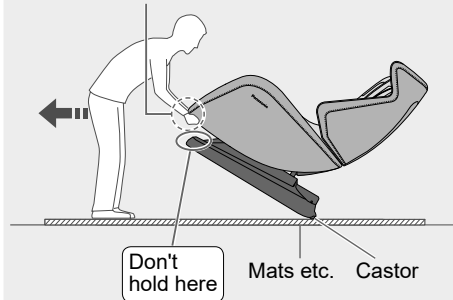
Before each massage (check the unit)

1 Check the surrounding area

- Make sure there are no objects, persons or pets around the unit.

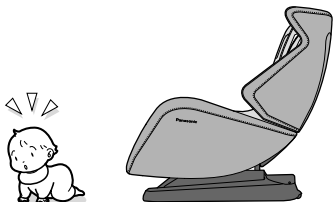
2 Flip the pillow

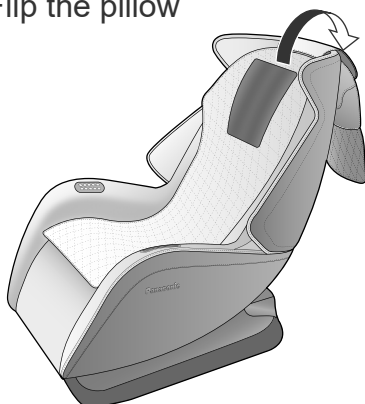
3 Make sure there are no foreign objects in the unit

If checking the leg massage section, turn the leg massage section down by hand, and check whether there are foreign objects in it.

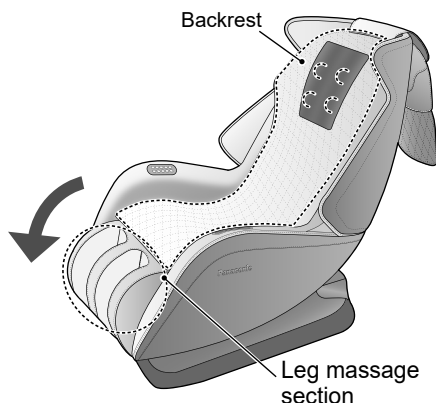
4 Check the power cord and power plug

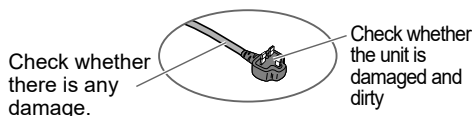⚠ WARNING

- Do not do anything that may damage the power cord or power plug. Do not scratch, modify, use near a heater, excessively bend, twist, pull, place heavy objects on top of, or bundle the cord or plug. Use damaged power cord or power plug may cause electric shock, short circuit or fire.
→ Contact the store where you bought the unit or the customer service center for repairs to the power cord or power plug.
- Dust the power plug on a regular basis. Remove dust using a dry cloth. (Otherwise it may lead to fire, due to damp shorting the insulation.)

5 Turn on the power switch.

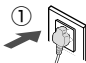- ① Insert the power plug into the power outlet. 

- ② Turn on the power switch.

- When the power plug is plugged in and the power switch on the back of the unit is turned on, the indicator light on the controller panel  starts flashing.

Settings


Before each massage (check the unit)

Before each massage (check the unit) (continued)

WARNING

- Do not insert or unplug the power plug with wet hands.
(Otherwise it may lead to electric shock.)
- Always use sockets and wiring devices at the correct rated value.
- Use 220 V 50 Hz alternating current.
(Cannot be used overseas, nor with a transformer.)
- Always insert the power plug fully into the socket.
(Otherwise it may lead to electric shock or fire due to overheating.)


- 6** Touch and confirm that the massage roller is in the original position

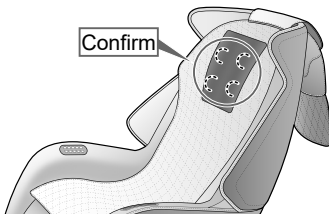When the massage rollers are not in the original position

- Press twice  on the controller panel to turn it back to the original position.

- 7** Make sure the fabric is undamaged

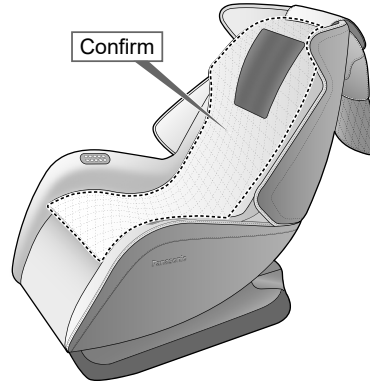WARNING


- Before use, be sure to check whether the covering fabric is damaged.
(If it is damaged, please stop using it immediately and entrust repair.)
Please do not use the unit when the fabric has been ripped, otherwise it may cause injury or electric shock.

- 8** Adjust the massage lounger as needed


When you want to use the leg massage section

Turn the leg massage section down and put your feet into the section

Turn the leg massage section down by hand

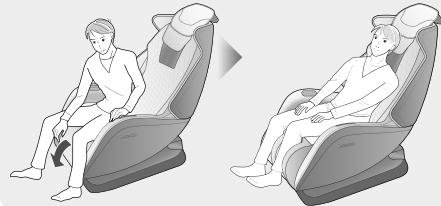Select menu and start massage

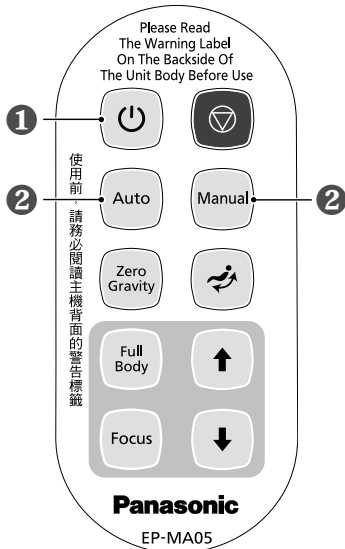1 Standby key

- After the power is turned on, the indicator light starts flashing (about 30 seconds) to wait for operation. After 30 seconds, the indicator light stops flashing (for not affecting customers to have a rest).

- Press for power-on
The indicator light is always on (standard time of the course is 15 minutes), and the backrest of massage lounge automatically tilts to the set angle (not the lowest position).

Note

15 minutes later, the course will stop, and the backrest of massage lounge will not reset automatically (for not affecting customers to have a rest).

- Press for shutdown
The indicator light starts flashing, and the backrest of massage lounge will reset to the chair state.

Note

15 minutes later, the course will stop, the switching course will not be re-timed.

2 Select massage course

Select a course from the "Auto" massage course and the "Manual" massage course and press it

Auto course


There are 3 courses to choose from, which can massage from shoulders to hips according to the preset massage course.

Manual course

Customers can set the massage techniques and massage parts according to their needs.

Before each massage (check the unit) \ Select menu and start massage

"Auto" massage course

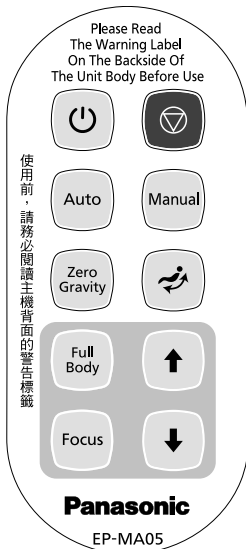Press  to enter the "Auto" massage course

 "Auto" massage course has three Pre-set course, which can be switched cyclically.

- "Pre-set course 1"
the red indicator light will be on. "Swedish" course.
- "Pre-set course 2"
the blue indicator light will be on. "Refresh" course.
- "Pre-set course 3"
the purple indicator light will be on. "Deep" course.

Course contents

Swedish

- When you want to eliminate the fatigue of the whole body, you should massage comfortably and gently from the shoulder to the hips by kneading and tapping, and squeeze and knead the calf softly and forcefully, so that your body can relax wholeheartedly.


Refresh

- When you want to relieve the stiffness of shoulder and waist, you should mainly massage and knead the back, and the massage rollers move at a slow speed, which can dredge the meridians and relax the body.

Deep

- In order to relieve the soreness and fatigue of the waist and hips the waist and hips can be better relaxed by tapping and rubbing, and the soreness can be thoroughly relieved by the squeezing and kneading of the calf air.

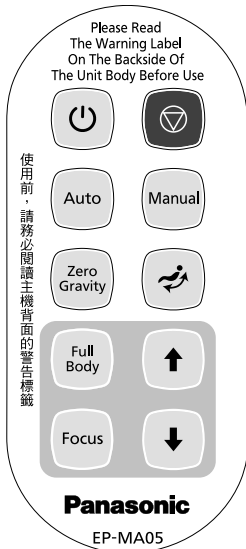"Manual" massage course

Press  to enter the "Manual" massage course

Select a course from "Full Body" or "Focus".


 "Full Body" massage course has three Pre-set course, which can be cycled by pressing , and can relax the whole body from shoulder to hip.

- "Pre-set course 1"
the red indicator light will be on. Mainly by kneading.
- "Pre-set course 2"
the blue indicator light will be on. Mainly by tapping.
- "Pre-set course 3"
the purple indicator light will be on. Mainly by kneading and tapping.

 "Focus" can select the part you want to massage, long press on  or  to adjust it up and down.

"Auto" massage course \ "Manual" massage course

Set up according to your favorite action

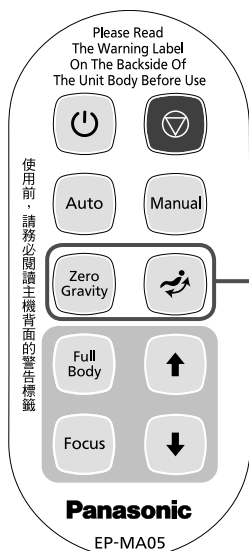Press  and the massage lounger will tilt automatically to the horizontal position.

- Incline the body to the best angle for "Zero Gravity" massage, so that the stress of the whole body can be balanced without any obvious stress points. At this time, other parts of the body can be better relaxed, which will be helpful to release the pressure and fatigue, and can also achieve better massage results.

Press  to enter constant "Swing" mode.

- Through natural and relaxed "Swing" massage, the body can be fully stretched, thus every part of the body can be massaged.

Music playing instructions

The massage lounger has built-in Bluetooth receiving module, which can be connected with Bluetooth mobile phones, tablet computers and other devices to play music, so that you can enjoy music while having massage. (EN21)

To ensure the security of Bluetooth connection, the password should be put in as follows.

Bluetooth ID: EP-MA05-XXXXXX

Bluetooth password: 2105

After completing the massage


1 Remove the feet from the leg massage section

For your safety, please put it away after the leg massage section is used.


Retract the leg massage section

- Please make sure there are no foreign material in the storage section.

Carefully remove the feet from the leg massage section



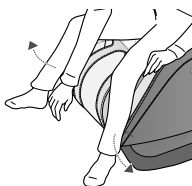
Push inward to the end with your hands

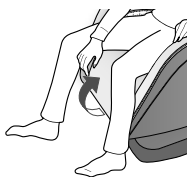WARNING

- This unit should not be used by children. Also, do not allow children to play on the unit. Do not allow children to climb on the leg massage section, seat, backrest or armrests. (Otherwise it may lead to accident or injury.)

2 Turn the power off

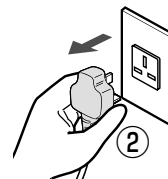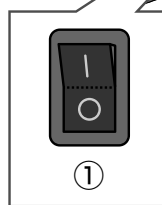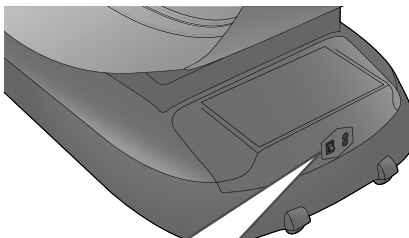- ① Press the power switch to the "O" off position
- ② Remove the power plug from the power outlet.

WARNING

- After each massage, push the power switch on the unit backside to the "O" (off) position. (To prevent breakage or injury from misuse by children.)

Timer function

After starting the massage, it will end automatically about 15 minutes later.

When you want to end it in the middle of the massage

Press 

- The massage rollers should be retracted and the backrest be raised.

When feeling anything abnormal

Press 

- Stop all actions immediately.
- Be careful when leaving the unit to avoid falling down.

Set up according to your favorite action After completing the massage

Cleaning and maintenance

- Never use chemicals such as thinner, benzene, alcohol, etc.

Synthetic leather/plastic parts

Normal maintenance


Use a soft dry cloth to wipe off dust and dirt.

If there is smear

- 1 Soak a soft cloth in water or a 3 % - 5 % solution of neutral kitchen detergent with warm water and squeeze well.

- 2 Wipe the surface.
(Clean and wipe cloth and synthetic leather)

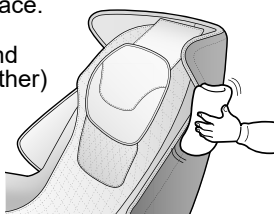- 3 Finally, rinse the cloth in water and wring it out, then wipe off all the detergent.

- 4 Wipe with a soft and dry cloth.
- 5 Let it air dry naturally.
(Do not rapidly dry off with a dryer, etc.)

If it is not easy to clean

Please try the above steps 1 by adopting the following maintenance methods.

<Synthetic leather>

Wipe with the available "melamine foam sponge" on the market with neutral kitchen detergent.


Note

- As synthetic leather component may be dyed by clothes, so please be careful when using the unit in jeans or colorful clothes.
- To avoid discoloration of synthetic leather, please put on a towel or other mat when using the massage lounge for hairdressing.


EN20

Fabric parts

- 1 Soak a soft cloth in water or a 3 % - 5 % solution of neutral kitchen detergent with warm water and squeeze well.
- 2 Use water or a neutral kitchen detergent and brush it on areas where the seat fabric is particularly soiled. Be careful not to brush the fabric too much to avoid damage. Detergent is not recommended for maintenance in the elastic fabric range. The scope of maintenance should be limited to the leather parts.)

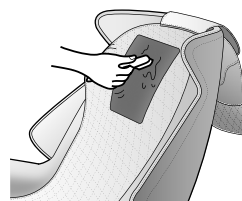- 3 Finally, rinse the cloth in water and wring it out, then wipe off all the detergent.

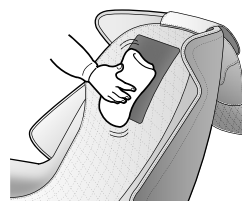- 4 Let it air dry naturally.

Note

- Since color may transfer from clothing to the seat fabric, be careful when using the unit with clothing such as jeans or colored clothes.

Q&A

Q Can persons paying frequent visits to hospitals because of diseases use the product?

A These persons should consult physicians before using the product.
Massage is a form of stimulation by tactile pressure in which pressure is applied to the muscles to relax them. According to different causes, massage may worsen the condition, so please consult with the attending doctor before using this product.

Q How tall of people can be suitable for using this product?

A Depending on the size or height, there will be some differences. Please use this product with general reference below.

- When users with a height of more than 170 cm use head massage, there may be insufficient contact of the head massage part or inability to massage the head, which could affect the massage effect.


For a taller person with knee bending when using the product:

Please choose either of the following method **A** or method **B**.

- A** Put down the backrest and lower the position of the leg massage section.
- B** Use the leg massage section to massage only the upper body.

Q Where can I see the manufacturing number of the product?

A It is below the rating sticker on back of the main body.

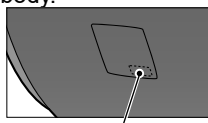Product serial number

Q A little concerned about the sound of the massage lounger when it works ...

A The sound produced when the massage lounger operates is usually caused by the moving structure.

Q Why does it stop after a period of time when using swing mode?

A Because the course is set to swing for 3 minutes and stop for 2 minutes, it is normal to stop for a period of time after swinging.

Q Does it matter that there is a difference in massage strength between the left and right massage rollers during massage?

A Due to the movement structure, it is possible to feel the difference in the strength between the left and right massage rollers during massage. An alternation tapping mechanism on the left and right massage rollers of the structure is used to simulate natural massage and create a more natural sensation. During the massage process, the specific positions of the left and right massage rollers in contact with the back may be slightly different, which could lead to the user's feeling that the massage strength on the left and right sides is different. This is not a malfunction.

Q What is the monthly power consumption when this product is used?

A According to the daily use of 30 minutes (15 minutes × 2 times), the electricity consumption of this product is about 1kW.h per month.

Q Are there any things needing special attention during moving?

A Pay attention to the following things during moving.

- For convenience during moving, tilt backrest of the product till the lowest position. During handling, do not sit on the product.

Turn off the power switch and pull out the power plug.

- If the original packing box is available, make sure the power line is not wound on the product, wind it and place it on the seat surface. When handling the product, keep the power line from falling off and rubbing with ground.
- If the original packing box is unavailable, properly keep other packing auxiliary materials and wrap the product with them to avoid scratch and vibration during handling.

Q How to connect and use Bluetooth audio of the massage lounger?

A After the massage lounger is turned on, turn on the mobile phone and open the Bluetooth to search for the massage lounger model [EP-MA05-XXXXXX] and connect it, and select music playing.

Troubleshooting


After discovering the fault, first check it and try the following methods.
If the fault still exists after trying, please contact the customer service center.

Problem	Cause and Remedy
<p>The massage rollers stop during the use of the unit</p>	<p>When the massage rollers encounter excessive external force during use, the massage rollers may stop moving for safety reasons. Please turn off the power switch, wait for about 10 seconds, and then turn on the power again. If the massage rollers stop again, please adjust your sitting position so that your back will not be completely attached to the backrest. (Attention should be paid to the relevant safety instructions, especially when used by people weighing 100 kg or above.)</p>
<p>The massage rollers do not come up to the shoulder or neck</p>	<p>If the massage ball does not reach the shoulder or neck position, you should adjust the button. Please refer to EN17 for specific operation instructions.</p>
<p>The height (strength) of the left and right massage rollers is different</p>	<p>An alternating tapping method is employed, which means that this phenomenon will naturally occur. There is nothing wrong with the unit.</p>
<p>The unit will not operate at all</p> <ul style="list-style-type: none"> ● Upper body (from neck to waist) massage ● Leg massage 	<ul style="list-style-type: none"> ● The power cord is disconnected. ● The power switch of the unit is not turned on. ● After pressing the  button on the controller, no course is selected.
<p>After massage, even if you press the  button again, you cannot re-enter the massage course</p>	<ul style="list-style-type: none"> ● It takes about 5 seconds for the course to stop after massage. Wait at least 5 seconds once all movement has stopped before pressing .
<p>The unit has been damaged</p>	<p>To prevent accidents, please contact the customer service center when necessary.</p>
<p>The power cord or power plug is abnormally hot</p>	

Sounds and sensations of the unit

The following action sounds and sensations will occur during use, but this is due to the structure of the massage lounger, nothing abnormal.

- The sound caused when the massage rollers are changed from "extended" to "retracted"
- The sounds and sensations when the massage rollers rub the wrinkles on the fabric during the movement
- The sound of rubbing, tapping, etc.
- Friction sound when massage rollers move
- Action sound when massage rollers go up and down
- Sound of the belt rotating
- Motor sounds
- The sound when the massage rollers change from "squeezing" to "stretching"
- Voice when seated
- Friction between armrest and pillow when backrest inclines
- Friction sound between massage roller and cloth
- Working sound of air pump during air massage operation
- Exhaust sound during air massage operation
- Working sound of solenoid valve under seat during air massage operation
- Shake of leg massage section
- Action sound of air
- When the backrest is tilted, the sound made by the friction between the pillow and the leather part on the back cover

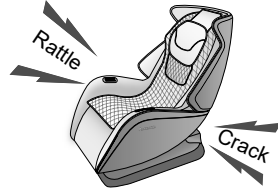When the sound becomes louder, that may suggest a fault/malfunction.


- At this time, please contact the customer service center for inspection and repair.

Specifications

Power supply		220 V ~ 50 Hz
Rated power		70 W
Massage range	Kneading width	122 mm - 47 mm
	Moving distance	86 cm
	Massage speed	Kneading width: 31±10 % times/min Tapping: 700±10 % times/min Moving distance: 3±10 % cm/s
Turn off automatically		Turn off automatically after approx. 15 minutes at most
Dimensions	When not reclining (length × width × height)	Approx. 117 cm × 68 cm × 106 cm
	When reclining (length × width × height)	Approx. 132 cm × 68 cm × 100 cm
Mass		Approx. 43 kg
Packing case dimension (D × W × H)		Approx. 146 cm × 76 cm × 96 cm
Weight (including packing case)		Approx. 53 kg
Upholstery		Synthetic leather
Bluetooth Version		Bluetooth 5.0

“The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Panasonic Corporation is under license. Other trademarks and trade names are those of their respective owners.”

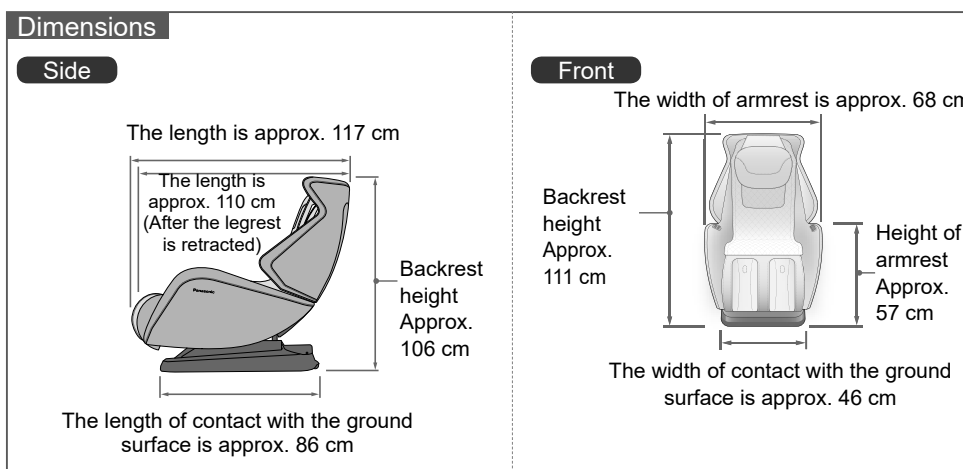For persons weighing 100 kg or more, the product operations may produce more noise and the fabric cover may wear faster.

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© Panasonic Corporation 2021