






Panasonic


Operating Instructions
(Household use) Sonic Vibration Toothbrush
Model No. EW-DA48

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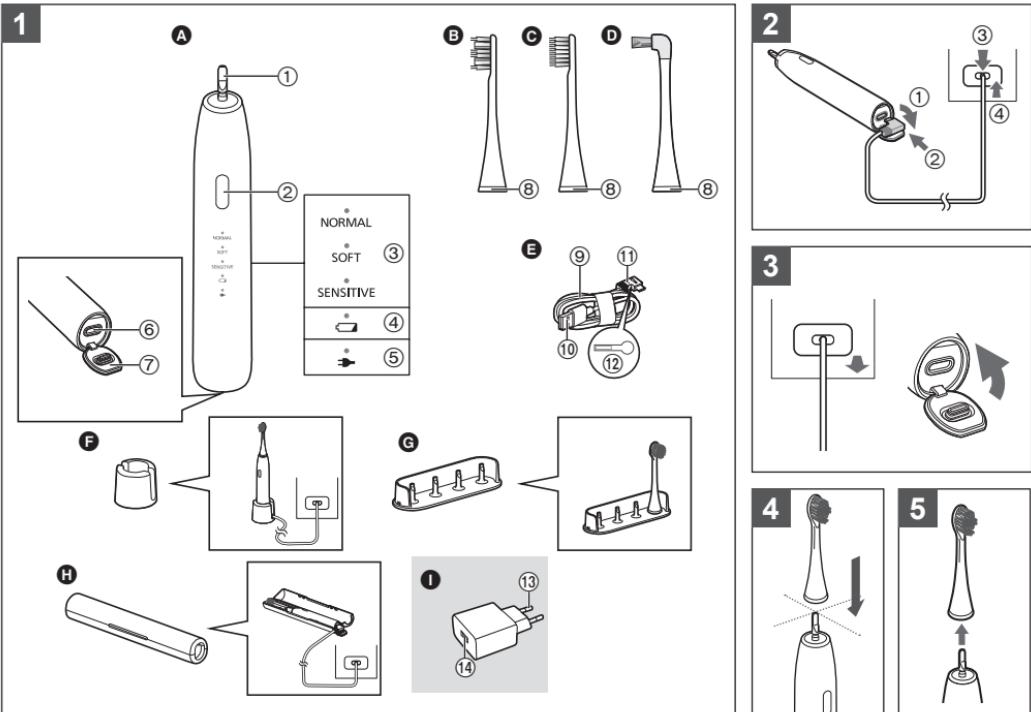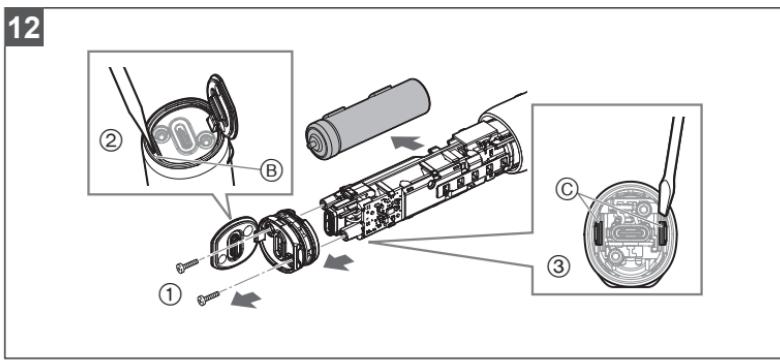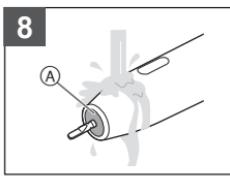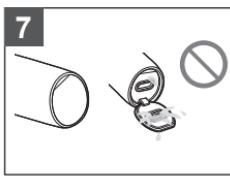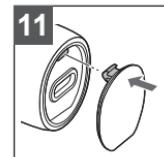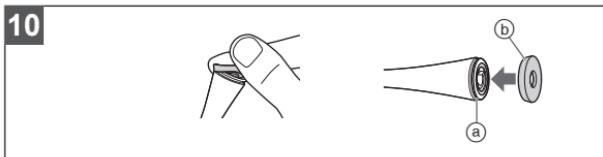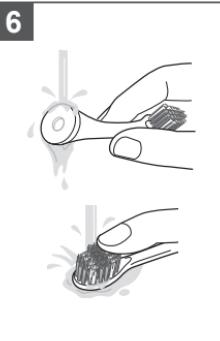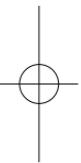Panasonic

English

Operating Instructions (Household use) Sonic Vibration Toothbrush Model No. EW-DA48

Contents

Safety precautions	6	Optional parts (Sold separately).....	15
Important information	9	FAQ	16
Parts identification	10	Troubleshooting	17
Before cleaning your teeth	10	Battery life	20
Features	10	Removing the built-in rechargeable battery	20
Charging	11	Specifications	20
Basic usage	12		
Tips on usage	13		
After use	14		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

The USB power adaptor is sold separately.

Safety precautions

Make sure to follow these instructions.

To prevent accidents, injuries or damage to property, please follow the instructions below.

■ The following chart indicates the degree of damage caused by improper operation.

 **DANGER**

Indicates a potential hazard that will result in serious injury or death.

 **WARNING**

Indicates a potential hazard that could result in serious injury or death.

 **CAUTION**

Indicates a potential hazard that could result in minor injury or property damage.

 **WARNING**

► **Power supply**

- Do not use in a way exceeding the rating of the household outlet or the wiring.
 - Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.

- Fully insert the power plug of the USB power adaptor securely.
- Failure to do so may cause fire or electric shock.

► **USB power adaptor and USB cable:**

- Do not immerse them in water or wash with water.
- Do not place them over or near water filled sink or bathtub.
- Never use the appliance if they are damaged or if the power plug fits loosely in a household outlet.
- Do not damage them.
 - Do not deface, modify, forcefully bend or pull them. Also, do not place anything heavy on, twist or pinch them.
 - Doing so may cause electric shock or fire due to a short circuit.
- Do not connect or disconnect them with a wet hand.
 - Doing so may cause electric shock or injury.
- Do not use them to charge other devices.

[Note that USB power adaptor (ES-RU1-02) can be used with compatible Panasonic products.]

 - Doing so may cause burns or fire due to a short circuit.
- Disconnect the USB plug when cleaning the main unit.
 - Failure to do so may cause electric shock or injury.
- Regularly clean them to prevent dust from accumulating.
 - Failure to do so may cause fire due to insulation failure caused by humidity.
 - Remove the USB power adaptor and wipe it with a dry cloth.

WARNING

- Always plug the USB cable to a power supply with a rated voltage of 5.0 V and a rated output current of 1.0 A which meets requirements of the Safety Extra Low Voltage (SELV) system.
 - Failure to do so may cause burns or fire due to a short circuit.
- Use the specified USB power adaptor and USB cable.
 - Failure to do so may cause burns, electric shock or fire due to a short circuit.

► In case of an abnormality or malfunction

- Immediately stop using the device and disconnect the USB power adaptor.
 - Failure to do so may cause fire, electric shock, or injury.
- <Abnormality or malfunction cases>
- The main unit, USB power adaptor or USB cable is deformed or abnormally hot.
- The main unit, USB power adaptor or USB cable smells of burning.
- There is an abnormal sound during use or charging of the main unit.
 - Immediately request inspection or repair at an authorized service center.

► This product

- This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, leave, or store in a high temperature environment.
 - Doing so may cause overheating, ignition, or explosion.

- Do not modify or repair.
 - Doing so may cause fire, electric shock, or injury. Contact an authorized service center for repair (battery change, etc.).
- Never disassemble except when disposing of the product.
 - Doing so may cause fire, electric shock, or injury.

► Preventing accidents

- Do not store within the reach of children or infants. Do not let them use it.
 - Doing so may cause an accident or injury due to accidental ingestion of removable parts.

► Note the following precautions

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
 - Failure to do so may cause an accident or injury.

CAUTION

► Preventing damage to teeth and gums

- Use the toothbrush with a general type of toothpaste only. Do not use specially formulated products such as ones enhanced for tar removal.
- Do not apply the brush head strongly against teeth or gums, or use it continuously against same spot.
- Do not apply the plastic part of the main unit or brush head to teeth.
- Do not let anyone who cannot operate the brush by themselves use it.
- People with low sensitivity in their mouth should not use the brush.
 - Doing so may cause injury or damage to teeth and gums.
- People who may have periodontal disease or who are concerned about other symptoms in their mouth should consult a dentist before use.
 - Failure to do so may cause injury or damage to teeth and gums.

► This product

- Never use for a purpose other than cleaning your teeth.
 - Doing so may cause injury.
- Do not share the brush head with your family or other people.
 - Doing so may result in infection or inflammation.

► Note the following precautions

- Do not drop or hit against other objects.
 - Doing so may cause injury.

- Do not allow foreign objects like metallic pins or lint on the USB power adaptor's power plug or USB plug.
 - Doing so may cause electric shock or fire due to a short circuit.
- Do not charge when the appliance socket is wet.
 - Doing so may cause burns or fire.
- Disconnect the USB power adaptor from the household outlet when not charging.
 - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

Handling of the removed battery when disposing

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- Never peel off the tube.
 - Doing so may cause overheating, ignition, or explosion.

WARNING

- After removing the rechargeable battery, do not keep it within reach of children and infants.
 - The battery will harm the body if accidentally swallowed.
If this happens, consult a physician immediately.

- If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.
 - Failure to do so may cause eye injury.
Consult a physician immediately.

CAUTION

- If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.
 - Failure to do so may cause inflammation.

Important information

- Charge the appliance correctly according to these Operating Instructions. (See page 11 "Charging".)
- This appliance contains a battery that is only replaceable by skilled persons. Please contact an authorized service center for the details of repair.

This symbol on the product means that not generated higher voltage than the safety extra low voltage (SELV). In which protection against electric shock relies on power supply.

Parts identification →P.2-1)

A Main unit

- ① Toothbrush shaft
- ② Power switch
(power off/on mode selection)
- ③ Mode indicator (blue)
- ④ Battery capacity lamp () (blue)
- ⑤ Charge status lamp () (red)
- ⑥ Appliance socket (Type-C)
- ⑦ Socket cover

Accessories

- F Main unit stand
- G Brush stand

Optional parts

- I USB power adaptor for general use (ES-RU1-02)
 - ⑬ Power plug
 - ⑭ USB port (Type-A)

B Total care brush

- Purple bristles
- ⑧ Identification ring

C Extra fine brush (compact)

- Blue bristles

D Extra fine point brush

E USB cable

- ⑨ Cord
- ⑩ USB plug (Type-A)
- ⑪ USB plug (Type-C)
- With overheat protection function
- ⑫ Overheat protection function mark

H Carrying case

Before cleaning your teeth

► This product

The Sonic Vibration Toothbrush is an electric toothbrush with sonic vibration* of approximately 31 000 brush strokes/minute.

* Sonic vibration in acoustic wave field

► Bleeding of the gums

When you first begin using the toothbrush, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the toothbrush for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the toothbrush and consult your dentist.

Features

► Overpressure protection sensor


To avoid damaging the teeth and the gums by brushing too hard, the toothbrush will temporarily stop its brushing action when the sensor detects excessive brush contact pressure, and alerts the user with a series of rapid beeps.

- This feature will not be available in the SENSITIVE mode.
- You can choose to turn the sensor function on or off. (it will be on by default)

How to turn the sensor on or off

Press and hold the power switch for 5 seconds while the brush is being charged.

- "Setting Change Successful" tones
 - When the sensor is turned from on → off: One long beep
 - When the sensor is turned from off → on: One short beep

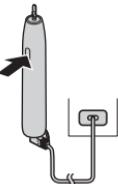► When to charge

The toothbrush will give lamp and tone signals to alert the user when the battery charge level gets low.

- When the toothbrush is switched on: the battery capacity lamp () will blink (once per second)
- When the toothbrush is switched off: A series of beeps will be given (Two beeps) (the lamp turns off after 5 seconds)

► Even when the toothbrush becomes completely discharged

- The toothbrush will be good for use for about 2 minutes after charging for 4 minutes. (Quick Charge feature)

(The available use time may vary depending on the ambient temperature and the manner in which the toothbrush is used.)

Charging (→P.2-2)

1 Open the socket cover and insert the USB plug firmly into the main unit and USB power adaptor.

- Wipe any water or toothpaste residue off the appliance socket area.

2 Insert the USB power adaptor firmly into a household outlet.

- Charging is completed after approx. 2 hours.

3 Disconnect the USB power adaptor from the household outlet when the charge status lamp () turns off. (→P.2-3)

- Close the socket cover.

Turns off when the toothbrush is fully charged. (2 hours)

► Battery runtime on a full charge

Mode	Battery runtime
NORMAL	
SOFT	Approx. 90 minutes
SENSITIVE	

(The available use time may vary depending on the ambient temperature and the manner in which the toothbrush is used.)

Notes

- The charging time may vary when using for the first time, or when it has not been in use for more than 6 months. Also, the charge status lamp () may not glow for a few minutes after charging starts, but it will glow if charging continues.
- The required charging time may vary depending on the starting charge level.
- Recommended ambient temperature for charging is 10 °C – 35 °C. When the ambient temperature is out of the recommended range, the battery may not charge properly or may take a longer time to be fully charged.
- If the charge status lamp () blinks quickly (twice per second), see page 17.


Basic usage

1 Attach the brush head to the main unit.

(→P.2-4)

2 Apply toothpaste.

There is a small gap between the brush head and the main unit (a few mm), but this is not a malfunction.

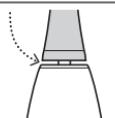
3 Put the brush head in your mouth and press the power switch to turn on the main unit.

12

► Select the desired mode

To switch the mode, press the power switch within 5 seconds.


- The device starts in the previous mode used when the power is turned on.

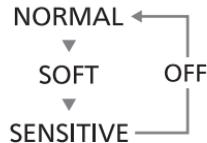
4 Lightly press and move the toothbrush against your teeth.

► Example brushing sequence

- front of the teeth
- biting plane of the teeth
- back of the teeth

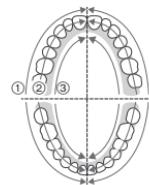The toothbrush will pause once at 30-second intervals to indicate the brushing time. It pauses twice at 2-minute intervals then starts again.

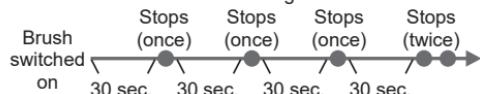5 Turn off the main unit before removing it from your mouth.

**Notes**

- The appliance may not operate in an ambient temperature approximately 5 °C or lower.
- Use the appliance in the appropriate mode.

Mode	Strength	Usage or desired care
NORMAL	Strong	To brush with standard strength level
SOFT		When you are using a powered toothbrush for the first time or want to brush gently
SENSITIVE	Weak	To brush very gently to avoid irritating the gums

- The toothbrush will temporarily stop its brush action when it senses excessive brush contact pressure against the teeth and the gums, and alerts the user with a series of rapid beeps. (See page 10.)

Tips on usage

Select the appropriate brush head in accordance with your objective and the area you are brushing.

For people with good, healthy gums**Total care brush**

This one brush head serves all your needs for gum care, plaque removal and stain control.

For people who have gum health concerns**Extra fine brush (compact)**

This brush head is specially designed to effectively remove deposit in narrow spaces around the teeth, such as periodontal pockets and the gaps between teeth.

To effectively clean hard-to-reach areas of teeth**Extra fine point brush**

This brush head is good for cleaning areas of the teeth that are hard to reach and are difficult to brush well, such as the back of the molars or where you have a crooked or misaligned tooth.

How to apply/move the toothbrush

Front surface and the biting surface of teeth

Periodontal pocket and the boundary between the teeth and gums

For the back of the front teeth or misaligned teeth

- Move toothbrush to match the angle of the teeth.


For the misaligned teeth

Back of the molars

Periodontal pocket/space between the tooth and the gum

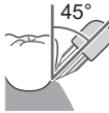After use

- Do not use thinner, benzine, alcohol, etc. to clean the appliance.

Doing so may cause failure, cracking or discoloration of the parts.

Remove the brush head. (→P.2-5)

Brush head

Wash the brush head under running water.

(→P.3-6)

- Do not wash with water hotter than 80 °C.
Doing so may damage the brush head.

Main unit

Close the cover and wash with water

Before washing the toothbrush, push the cover down to see it is tightly closed. (→P.3-7)

Wash the main unit under running water.

(→P.3-8)

- Do not damage or rub the waterproof rubber seal (Ⓐ) with hard materials.
Do not leave the main unit under water.
- Do not wash with water hotter than 40 °C.
Doing so may damage the main unit.

Wipe off water from the brush head and main unit for storage. (This is to keep the toothbrush clean and to reduce malfunction risks.)
• Do not charge the toothbrush when it is wet.

Main unit stand/brush stand/ carrying case

Wash and wipe off water. (→P.3-9)

Replacing the identification rings (→P.3-10)

► Removing

Slide up with your thumb.

► Attaching

Insert by covering the indentation (Ⓐ) in the brush head with the inside of the identification ring (Ⓑ).

If the socket cover comes off

Push the cover onto the appliance socket straight down. (→P.3-11)

Optional parts (Sold separately)

- Brushes are consumables.
- The brush should be replaced with a new one every 3 months.

Part number	Replacement Brush for EW-DA48
WEW0810	Total care brush
WEW0800	Extra fine brush (compact)
WEW0861	Extra fine point brush

- Only one of each brush head is included in the main unit's package.

ES-RU1-02 | USB power adaptor

FAQ

English

Question		Answer
While using	My toothbrush temporarily stopped with rapid beeps. Why?	<p>The overpressure protection sensor has been activated. Open your mouth lightly and apply the brush to your teeth with gentle pressure.</p> <p>If the beeps bother you, you can choose to turn the sensor off. (See page 11.)</p>
Battery	Why does my toothbrush temporarily stop or reduce the brush action about every 30 seconds?	<p>The toothbrush will notify you of the elapsed brushing time every 30 seconds.</p> <p>(See page 12 for description of the brushing time notification function.)</p>
	Will the battery deteriorate if not used for an extended period of time?	<p>The battery will deteriorate if not used for over 6 months, so make sure to perform a full charge at least once every 6 months.</p>
	Can I charge the appliance after every use?	<p>The battery used is a lithium-ion, so charging after every use will have no effect on battery life.</p>
While charging	Can I use a commercial USB cable or USB power adaptor to charge my toothbrush?	<p>Although the toothbrush may be rechargeable with some products other than our designated accessories, please use our designated charging accessories listed below as they are fully verified to be compatible with the toothbrush.</p> <ul style="list-style-type: none">• Included USB cable: Has overheat protection feature• Separately available USB power adaptor: Splash-resistant design (not fully waterproof; avoid wetting the toothbrush.)

Troubleshooting

Problem	Possible cause	Action	
While charging	<p>The charge status lamp () does not glow when charging starts.</p> <p>The charge status lamp () is blinking rapidly. (Twice every second)</p>	<p>► You have just purchased the toothbrush or have not used it for a while.</p> <p>► The USB cable is not correctly connected with the main unit and/or the USB power adaptor.</p> <p>► The USB power adaptor is not correctly connected to the household outlet.</p> <p>► The temperature is too low or high while charging.</p>	<p>► The charge status lamp () may not glow for several minutes after starting to charge it, but it will glow if you continue to charge it. If the charge status lamp () does not glow after a while, consult with an authorized service center for repair.</p> <p>► Make all connections correctly and check that the charge status lamp () comes on. (See page 11.)</p> <p>► Charge the toothbrush within the recommended charging temperature of 10 °C to 35 °C.</p> <p>► If the charge status lamp () blinks even when charging within the recommended temperature range, please contact an authorized service center for repair.</p>

Problem	Possible cause	Action
While charging Toothbrush does not charge. The toothbrush can only be used for a few minutes, even if it is fully charged.	<ul style="list-style-type: none">▶ The USB cable is not correctly connected with the main unit and/or the USB power adaptor. The USB power adaptor is not correctly connected to the household outlet.▶ The wet appliance socket activated the overheat protection feature of the USB cable to suspend the charging action.▶ The rechargeable battery has reached its operating life. (Approx. 3 years)	<ul style="list-style-type: none">▶ Make all connections correctly and check that the charge status lamp () comes on. (See page 11.)▶ Allow the appliance socket to dry before recharging.▶ Replace the battery.▶ Contact the authorized service center where purchased.

If the problems still cannot be solved, contact the store where you purchased the unit or a service center authorized by Panasonic for repair.

Battery life

English

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge.
(The operating life of the battery may differ depending on usage or storage conditions.)

Removing the built-in rechargeable battery (→P.3- 12)

Remove the built-in rechargeable battery before disposing the toothbrush.

The battery is to be disposed of safely.

Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing the toothbrush, and must not be used to repair it. If you dismantle the toothbrush yourself, it will no longer be waterproof, which may cause it to malfunction.

- Press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Follow the procedure below, using a screwdriver to disassemble the appliance.

Perform steps ① to ③ and lift the battery, and then remove it.

② Insert the tip of a flat screwdriver into the groove (Ⓐ) in the bottom cap, and pry the cap open with the screwdriver.

③ While pushing at the toothbrush shaft (→P.2-②-①), push the two battery holder tabs (Ⓑ) inward with the tip of the screwdriver to disengage the tabs. Afterward, remove the battery.

- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

For environmental protection and recycling of materials

This toothbrush contains a Li-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Output rating of applicable adaptor (sold separately)	5.0 V--- 1.0 A In compliance with the Safety Extra Low Voltage (SELV) system
Charging time	Approx. 2 hours (Room temperature: 10 °C – 35 °C)

This product is intended for household use only.

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© Panasonic Corporation 2025
Panasonic Corporation